

MICRO MOTORS **INTRAFLOW ANESTHESIA DELIVERY SYSTEM**



Introducing a new era in patient comfort and dentist productivity

The IntraFlow™ Anesthesia Delivery System is the latest technology in pain management for the dental professional, and a major breakthrough in fear management for the anxious dental patient. The IntraFlow System uses patented intraosseous injection technology in a simple, one-step technique that provides immediate pain relief and profound pulpal anesthesia in about one minute.

PATIENT COMFORT

- No lingering numbness to the cheek, lips or tongue
- Patient can immediately resume normal activities, talking, eating, drinking
- Facilitates more flexible scheduling
- Relieves patient anxiety regarding "The Shot"

DENTIST PRODUCTIVITY


- One step, easy to use
- Saves time, no waiting for patient to get numb
- Profound and reliable anesthesia
- Work in multiple quadrants during a single visit
- Reduces the need for multiple follow-up visits

Pro-Dex
Micro Motors™

"Extraordinary Innovation, Extraordinary Results"

For more information please call Pro-Dex, Micro Motors at: 800-562-6204 or contact your authorized distributor.

SPECIFICATIONS FOR INTRAFLOW ANESTHESIA DELIVERY SYSTEM

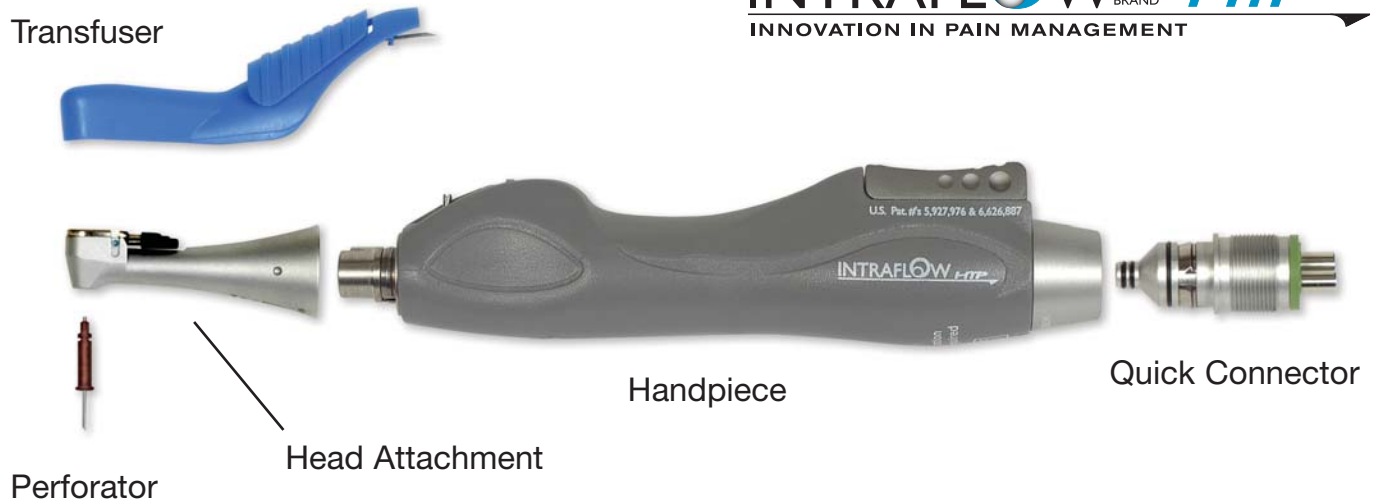

The IntraFlow™ Brand Anesthesia Delivery System is a product of Pro-Dex, Inc. It is a Class I medical device and regulated under Federal Code of Regulations Title 21, Parts 800, 801, 803, 806, 807, 820, 821, 860 and 872. Refer to the IntraFlow™ System Instruction Guide for a complete review of recommended uses and contraindications. Protected under U.S. Patent Laws.

U.S. Patent #'s : 5,927,976 – Medical injection device and method & 6,626,887 – Hard tissue drug delivery device and method.



Authorized EC Representative:
European Compliance Services, Ltd.
Oakdene House, Oak Road
Watchfield, Swindon
Wilts, SN68TD, UK

HANDPIECE & HEAD ATTACHMENT

Length	• 7.75" / 19.69 cm
Maximum Width	• 1.2" / 3.05 cm
Weight	• 4.8 oz / 135.4 g
Power	• Compressed air / 35 PSI (dynamic)
Adapter Type	• 4-Hole Midwest
Motor	• Variable speed drive
Anesthetic	• Standard (1.8ml cartridge)

PERFORATION & INJECTION SUBSYSTEM

Perforator	• 24 gauge (0.020" diameter) • Stainless Steel beveled needle
Transfuser	• ABS Shell & Slider • 20 gauge Stainless Steel Cannula

ORDERING INFORMATION

IF-SKT-004	IntraFlow™ System Introduction Kit
IF-HPC-001	IntraFlow™ System Handpiece
IF-ATT-003	IntraFlow™ System Head Attachment
IF-TRF-007	IntraFlow™ System Transfuser Assemblies - 20 count
IF-TRF-008	IntraFlow™ System Transfuser Assemblies - 80 count
IF-QDC-001	IntraFlow™ Quick Disconnect
IF-ACC-002	IntraFlow™ Perforator Recapping Stand
IF-ACC-003	IntraFlow™ Head Attachment Soaking Tub
IF-SIG-001	IntraFlow™ System Instruction Guide
IF-SCD-001	IntraFlow™ System CD-ROM

Note: Specifications are subject to change without notice.



151 E. Columbine, Santa Ana, CA 92707
t: 800-562-6204 f: 714-546-1109
www.pro-dex.com

**"Extraordinary Innovation,
Extraordinary Results"**